BS 38

THERMOMETER-SET

thermometer set • set de thermomètres • set termometro • kit de termómetros

Service-Hotline: 0180 5 007514

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GB USER MANUAL

Thank you for purchasing the H+H digital thermometer set.

Please retain the receipt with the original date of purchase for your records. In the case of service and warranty issues, you will always be requested to submit a copy of this receipt. Please contact your H+H dealer regarding warranty claims.

The user manual provides you with all information to ensure safe, easy use of the thermometers. Before using the thermometer set, we request you read this manual thoroughly. Please pay special attention to the safety instructions.

CONTENTS OF PACKAGE

- Digital flex-tip thermometer
- Digital pacifier thermometer
- Digital bath and room thermometer

SAFETY INSTRUCTIONS

- Only use the flex-tip and pacifier thermometer to measure human body temperature.
- Before each use, check the teat of the pacifier for signs of wear or cracks and pull it to ensure that it is firmly attached to the shield. Do not use it if a defect is found or if it is loose.
- Ensure an adult supervises the use of the thermometers.
- Store the thermometers out of reach of children, they are not toys.
- Temperature measurement is no substitute for consulting your physician.


CARING FOR THE THERMOMETERS

- Do not store in direct sunlight or at high temperatures.
- Do not disassemble the thermometers.
- Never use boiling water or a chemical solvent to clean the thermometers.
- Do not wash in dishwasher.
- Do not put the thermometers in a microwave or microwave-sterilizer.

BEFORE THE FIRST USAGE

Clean the digital flex-tip thermometer and the pacifier thermometer before first use.

NOTES ON DISPOSAL

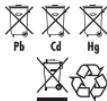Old units, marked with the symbol as illustrated, may not be disposed of in the household rubbish. You must take them to a collection point for old units (enquire at your local authority) or the dealer from whom you bought them. These agencies will ensure environmentally friendly disposal.

THE FLEX-TIP THERMOMETER

The digital flex-tip thermometer can be used to measure your baby's temperature rectally (rectum) or axillary (armpit).

USING THE FLEX-TIP THERMOMETER

1. Remove the protection cap from the thermometer.
2. Press the On/Off button to switch the thermometer on. „188.8E“ appears in the display, confirming that the thermometer is working.
3. Choose how the temperature should be measured. Rectal measurement is the most accurate and suitable for both infants and children. are not recommended for accurate measurements.

- When used rectally, the temperature sensor must be inserted 1 to 2 cm.
Vaseline can be used for lubrication.
The measuring time is 30 seconds.
The normal body temperature is about 37°C - 98.6°F.
- When used axillary (armpit), insert the temperature sensor in the armpit.
Press and fix with the arm.
The measuring time is 3-5 minutes.
The normal body temperature is about 36.5°C - 98.1°F.

i Note: The body temperature is different for all humans and changes during the day. You are therefore advised to measure the temperature at a specific time in the morning.

4. Start the measurement when the screen displays „Lo °C“ or „Lo °F“. During measurement, the symbol „°C“ or „°F“ will flash while the temperature increases. When a maximum value is reached, the symbol „°C“ or „°F“ ceases to flash. A buzzing sound that lasts about 10 seconds indicates that the measurement is completed.

10 Acoustic signals at 10 second intervals means that the measured temperature is normal (below 37.8 °C).
3 sets of short beeps at 10 second intervals means a possible fever (measured temperature above 37.8°C).

- When the display reads „Lo °C“ or „Lo °F“, the measured temperature is below 32.0°C/90°F.
 - When the display reads „Hi °C“ or „Hi °F“, the measured temperature is above 43.0°C/109.9°F
5. Press the On/Off button to turn the thermometer off.

i Note: The thermometer switches off automatically after 10 minutes. Switch the thermometer off directly after use to prolong battery life.

RECALLING THE LAST TEMPERATURE

Press and hold the On/Off button for about 2 seconds when switching the thermometer on. The last temperature measured will be displayed. This is indicated by the letter „M“ displayed under the degree symbol. After 3 seconds, the display reads 37.0°C and is ready to take a new measurement.

CLEANING

The thermometer is waterproof and can be cleaned with water, sterilization lotion or a soft cloth.

REPLACING THE BATTERY

The battery must be replaced when the display becomes indistinct or  appears.

To replace the battery:

1. Open the battery cover, the cover at the top of the thermometer.
2. Remove the old battery.
3. Insert a new battery and replace the battery cover.

SPECIFICATIONS

Measuring range:	32.0°C – 42.9°C (89.6°F – 109.2°F)
Measuring time:	approx. 3 – 5 minute
Accuracy:	± 0.1°C (±0.2°C) for 35.0 – 42.0°C (95°F – 107.6°F) and ±0.2°C (±0.4°F) for remaining range at an ambient temperature of 22°C (71.6°F)
Display:	LCD (display measuring unit: 0.1°C / 0.1°F)
Battery:	1.5VDC (LR41) alkaline
Power	consumption: 0.15 mW
Dimensions:	129mm(L) x 21mm(W) x 12.5mm(H)
Battery life:	approx. 250 hours or 1000 applications
Weight:	approx. 12 g. (with battery)

THE PACIFIER THERMOMETER

The most convenient way to measure your baby's temperature is with the thermometer designed as a pacifier.

Note: Oral (mouth) measurements are often about 0.5 °C or 1.0 °F lower than body temperature readings (rectal measurement).

USING THE PACIFIER THERMOMETER

1. Remove the protection cap from the pacifier.
2. Press the On/Off button. „188.8E“ appears in the display, confirming that the thermometer is working.
3. Insert the pacifier in the baby's mouth.

The measurement starts when the screen displays „Lo °C“ or „Lo °F“. During measurement, the symbol „°C“ or „°F“ will flash while the temperature increases. When a maximum value is reached (after approx. 3 minutes), the symbol „°C“ or „°F“ ceases to flash. A buzzing sound indicates that the measurement is finished.

Note: It is recommended to leave the pacifier thermometer in the baby's mouth for 6-8 minutes to allow the pacifier to warm up and to obtain a more accurate measurement.

- When the display reads „Lo °C“ or „Lo °F“, the measured temperature is below 32.0°C/ 90 °F
- When the display reads „Hi °C“ or „Hi °F“, the measured temperature is above 43.0°C/ 109.9 °F

Note: Drinking or eating prior to use may affect the temperature reading.
Wait half an hour before taking a measurement.

4. Press the On/Off button to switch the pacifier thermometer off.

Note: The thermometer switches off automatically after 10 minutes.
Switch the thermometer off directly after use to prolong battery life.

RECALL THE LAST TEMPERATURE

To recall the last temperature:

Press and hold the On/Off button for at least 2 seconds while switching the pacifier thermometer on. The last temperature measured appears in the display. This is indicated by the letter „M“ displayed under the degree sign.

CLEANING

Wash the pacifier thermometer after use with luke warm water and a mild soap.

BATTERY

When the battery is almost exhausted, the symbol appears on the display. This means that it is time to dispose of the unit.

Note: Do not use the pacifier thermometer as a normal pacifier when the battery is empty.

SPECIFICATIONS

Measuring range:	32.0°C – 43.9°C (89.6°F – 109.9°F)
Accuracy:	± 0.1°C (±0.2°F) for 35.0 – 38.0°C (95°F – 100.4°F) ambient temperature of 22°C (71.6°F)
Display:	LCD (display unit of measurement: 0.1°C / 0.1°F)
Power consumption:	0.15 mW
Dimensions:	40mm(L) x 43mm(W) x 53mm(H)
Battery life:	about 200 hours in continuous operation (approx. 2 years if used 10 minutes a day)
Weight:	approx. 14 g. (with battery)

THE BATH AND ROOM THERMOMETER

Use the bath and room thermometer to conveniently measure the temperature of the baby's room or the bathwater.

Note: Your baby will feel most comfortable in bath if the water temperature is between 36.5 °C and 38 °C. A temperature of 39 °C and above is too hot and your baby could be burned! Babies feel most comfortable when sleeping in rooms with a temperature of about 18 °C .

USING THE BATH AND ROOM THERMOMETER

Just let the thermometer float in the bathwater or let it stand on a shelf in the baby's room to start measuring.

Note: Mix the water thoroughly before measuring the water temperature or put the thermometer in the tub while filling it. This will ensure a fast and accurate reading of the water temperature.

- The display will read „COLD“ when the measured temperature is below 33 °C. Bathwater feels most comfortable at a temperature between 36.5 °C and 38 °C.
- Warm water alert: the red warning LED will start to flash when the measured temperature is 39 °C or above. The display will then flash „HOT.“

REPLACING THE BATTERY

The battery power starts to weaken after about 2 years. You can then replace it. To replace the battery:

1. Push the measurement unit out of the rubber housing
2. Remove the battery cover located at the back of the measurement unit using a coin or something similar.
3. Carefully remove the empty battery with a pair of tweezers.
4. Insert the new battery and replace the battery cover.

SPECIFICATIONS

Measuring range:	0°C – 50°C (32°F – 122°F)
Accuracy:	± 1°C (2°F)
Display:	LCD (display unit of measurement: 0.1°C / 0.1°F)
Battery:	1x CR2032
Dimensions:	40mm(L) x 43mm(W) x 53mm(H)
Battery life:	approx. 1 year (based on 1 hour of use a day)

GUARANTEE

In the case of a defect, please return the device together with the purchase receipt and original packing material to the point-of-sale.

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Olympia Business Systems Vertriebs GmbH