- DE** **AT** **CH** **Bedienungsanleitung**
- GB** **IE** **Operating Manual**
- FR** **Mode d'emploi**
- IT** **Istruzioni d'uso**
- ES** **Manual de instrucciones**
- NL** **Gebruikshandleiding**

BS 20

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Operating Manual

Please read and observe the following information and keep this operating manual safe for future reference!

Safety Instructions

⚠ WARNING!

- Risk of suffocation through small parts, packaging and protective foil!
- ➔ Keep the product and packaging out of reach of children!

⚠ CAUTION!

- Risks to health and the environment from batteries!
- ➔ Never open, damage or swallow batteries or allow them to pollute the environment. They may contain toxic, ecologically hazardous heavy metals.
- Risks to health from batteries!
- ➔ Keep batteries out of reach of children. Do not swallow batteries. If a round cell or battery is swallowed, contact medical assistance immediately.
- Property damage through overloading!
- ➔ Do not weigh any objects which are heavier than the maximum load stipulated in the technical data. Doing so could damage the scales permanently.


Intended Use

The digital baby scales, subsequently referred to as scales, serve to determine weights up to max. 20 kg. Any other use is considered unintended use. Unauthorised modification or reconstruction is not permitted. Under no circumstances open the scales or complete any repair work yourself.

Package Contents

- Baby scales
- Operating manual
- 2 batteries (AAA)

Operating Elements

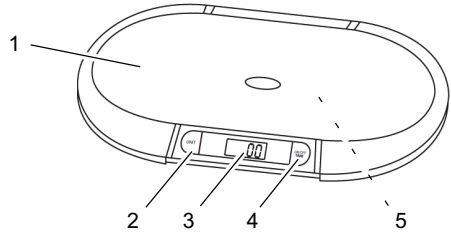- 1 Contact surface
- 2 Button for unit of weight
- 3 Display
- 4 Buttons for
 - a) Switching on and off,
 - b) Taring
- 5 Battery compartment

Inserting the Batteries

The scales are operated by means of two AAA batteries.

- ➔ Open the battery compartment cover on the underside of the scales.
- ➔ Insert the batteries (supplied with the device) in the battery compartment [5], ensuring correct polarity.
- ➔ Close the battery compartment cover on the underside of the scales.

i When the batteries are almost empty, the message **Lo** appears in the display. Change the batteries shortly after.

Operation

Switching on

- ➔ Set the scales down on a horizontal, stable base.
- ➔ Press the button [4] for 2 seconds. ➔ The scales are switched on.

Setting the unit of weight

The scales can display the weight in **kg** (kilogramme), **lb** (pound) or **st** (stone).

- ➔ Press the button [2]. ➔ The display [3] switches between the units of weight.

Taring the scales

If you weigh a baby or infant (or another object) using an aid (such as a towel), the scales must be tared.

i Taring means that the weight of the aid is taken into account for the weighing

process and the scales with the aid on it are set to zero, which ensures only the net weight of the baby or infant is determined.

- ➔ Place or lay the aid centred on the contact surface [1].
 - ➔ The display [3] indicates the weight of the aid.
- ➔ Press the button [4]. ➔ The display [3] then reads the value **0.00**.

Weighing

- ➔ Place the object to be weighed (with the aid previously tared, if necessary) centred on the contact surface [1]. ➔ The display [3] reads the weight measured.
- i** After a short time during which the weight does not change, the scales switch off automatically to save the batteries.
- i** If the object weighed is heavier than 20 kg, the message **O-Ld** appears in the display [3]. The display indicates the scales are overloaded.

Technical Data

Batteries	- 2 x AAA batteries
Load (max.)	- 20 kg
Weight	- 1.6 kg
Dimensions (H x W x D)	- 550 x 320 x 30 mm
Weighing accuracy	- 10 g

Care Instructions

- ➔ Clean the housing surfaces with a soft, lint-free cloth.
- ➔ Do not use and solvents or cleaning agents.

Disposal

In order to dispose of your device, take it to a collection point provided by your local public waste authorities (e.g. recycling centre). According to laws on the disposal of electronic and electrical devices, owners are obliged to dispose of old **electronic and electrical devices** in a separate waste container. The symbol indicates that the device must not be disposed of in normal domestic waste!

Risks to health and the environment from batteries!

Never open, damage or swallow batteries or allow them to pollute the environment. They may contain toxic, ecologically hazardous heavy metals. You are legally obliged to dispose of batteries and power packs at a battery retailer or in the corresponding containers provided at collection points by local public waste authorities. Disposal is free of charge. The symbol indicates that the batteries must not be disposed of in normal domestic waste and that they must be brought to collection points provided by local public waste authorities.

Packaging materials must be disposed of according to local regulations.

Guarantee

Dear Customer,
we are pleased that you have chosen this equipment.
In the case of a defect, please return the device together with the receipt and original packing material to the point-of-sale.

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Olympia Business Systems Vertriebs GmbH